



SERENITY

ZA-151

Massage Chair

User 's Manual



www.ZenAwakening.com

GALAXY HOME DESIGN
TULSA, OK
918-794-1539



www.ZenAwakening.com

Please read SAFETY INSTRUCTIONS carefully before using the chair.

SERENITY ZA-151 FUNCTION INTRODUCTION

1. With China national patent RK02 massage machine, the model can provide different massage combinations of over ten massage techniques such as kneading, rolling, tapping, shiatsu and knocking on head, neck, shoulder, waist and buttocks by the brand new mechanical programs to stimulate points, eliminate fatigue and promote blood circulation.
2. Air pressure massage on shoulder, arm, seat, leg and foot..
3. Leg-rest design which can extend or retract within 20cm, which is suitable for people in different height.
4. Brand new simple appearance design.
5. Back-rest and leg-rest is controlled by actuator to adjust their angles.
6. New human massage techniques and programs.

- ✦ 1. Thanks for choosing this SERENITY ZA-151 massage chair.
- ✦ 2. Please read this manual carefully before using the massage chair.
- ✦ 3. Please keep this manual in a safe place for further reference.
- ✦ 4. Please pay attention to SAFETY INSTRUCTIONS.

✦ Note: We reserve the right for design modification. It is subject to change without notice. The pictures in this manual are only for reference. Please refer to the actual product.

Content

Preparation	Safety instructions 1
	External structure parts 5
	Internal structure parts 6
	Mechanical components diagram 7
	Mechanical components list 8
Installation	Installation 9
Massage Operation	The adjustment before massage 13
	Manual control usage 15
	Operation instruction & function illustration 17
	After massage 19
Maintenance	Clean and maintenance 20
Others	Troubleshooting guide 21
	Technical data 22

Please keep this manual in a safe place for future reference.

- ★ Please read this manual carefully to make sure the installation and operation of the massage chair is correct.
- ★ “Safety Instructions” is provided in both Warning and Caution for important safety information.

Both Warning and Caution are important safety caution, please pay more attention.

 Warning	This sign indicates warning of causing possible death or hurt to body.
 Caution	This sign indicates caution of causing possible hurt or damage to things.

Symbol meaning

	This sign indicates forbidden behavior.		This sign indicates self disassembly is forbidden.
	This sign indicates must follow instruction to operate.		This sign indicates correct operation of cable.

- ★ The following should be kept for any further references.



Warning

	<ul style="list-style-type: none"> ● Osteoporosis patient; ● People with heart disease including embedded pacemaker or other medical electronic device inside the body; ● For pregnant women or people not healthy, please consult doctor before use; ● Kids and the drunk people are not allowed to use the chair; ● Please don't sleep in the chair. 	People who are not suitable to use the chair
	<ul style="list-style-type: none"> ● Make sure there are no kids or pets or other objects around (under, behind or in front of) the chair at the moment of backrest, leg-rest up/down or leg-rest retraction/extension from reclining position. ● During operation, if you have body abnormal, please stop the operation and go to consult your physician. ● Please do not take strong massage on the head. ● Do not massage belly or knees, be careful when massage neck. ● Children or pets are not allowed to play on the chair especially during massage. Standing or sitting on the back or armrest is not allowed either. ● Those who move around difficultly or can't express themselves must be taken care of strictly and continuously during massage. 	Accident/injury
	<ul style="list-style-type: none"> ● Insert the plug of the chair into correct and suitable power socket. Make sure the plug is completely inserted into the socket to reduce short circuit or fire risk. ● Unplug the power when not in use, stop all massage functions and turn off on/off key before unplugging. 	Fire and electric shock
	<ul style="list-style-type: none"> ● Massage cushions of the chair are not allowed to fold for storage. ● The machine is forbidden for no-load running for a long time. 	Accident
	<ul style="list-style-type: none"> ● Power cord, plug, socket can't be used in broken or loose contact condition. ● Don't use other electric power except rated voltage in case of fire, electric shock or any breakdown. ● Any broken, re-process, knotted, twisted power cord or pressed under heavy cargo is not allowed to use. 	Fire and electric shock
	<ul style="list-style-type: none"> ● If the chair works dysfunctional, or is broken, splashed with water or falls into water, please have it repaired by professional repairman or send the chair to assigned service agent. Self-modification, disassembly or maintenance is forbidden. ● If power cord or plug is damaged, please don't replace it yourself, please contact us or appointed dealer for replacement. ● To avoid electric shock or body hurt, please don't open or disassemble any part of the chair, there is not any attached part of the chair for maintenance by customers. ● Check leather or cloth material on the back and other areas is broken or not. If broken, even though it is small, please stop using the chair, unplug the power and have it repaired 	Disassembly and maintenance

**Caution**

	<ul style="list-style-type: none"> ● Please operate the chair in a flat or level place 	To avoid sudden falling.
	<ul style="list-style-type: none"> ● Position the chair 40cm from the wall to allow for sufficient space to recline up and down. 	To make sure the back lying.
	<ul style="list-style-type: none"> ● The maximum suggested operating time is 20mins. 5mins for each position beginning with the slight- step massage. Along continuous massage in the same position may cause uncomfortable feeling. 	Too long time massage is bad for your body.
	<ul style="list-style-type: none"> ● Before sitting on the chair to avoid injury please inspect everything to confirm that the massage heads are in the correct position. 	To avoid accidental injury.
	<ul style="list-style-type: none"> ● Please make sure that the power switch is at the "OFF" position before inserting the plug into the socket. 	
	<ul style="list-style-type: none"> ● Turn off all controls before removing plug from outlet. 	To avoid the accident when turn on the power.
	<ul style="list-style-type: none"> ● Make sure the plug completely insert into the socket. 	To reduce short circuits and any fire hazard.
	<ul style="list-style-type: none"> ● Checking whether any damage on wire, make sure no water or rubbish on wire. 	To avoid electric shock in any leakage or breakdown.
<ul style="list-style-type: none"> ● Installing earth wire in correct way. The plug must be inserted into socket with earth connect, the earth connect couldn't be bend or cut off. 		
	<ul style="list-style-type: none"> ● Don't press with hands or step on the massage heads with feet. ● Don't put hands between the working massage heads. ● Don't put hands or feet in the gap of mechanical parts. ● Don't drop or embed anything into the slot of massage chair. ● Don't stand or sit or put anything on the backrest, armrest or leg-rest. ● Don't stand on the chair during operation. ● Don't move the chair during operation. ● Don't pull out the plug or turn off the power during operation. ● Don't put your hands or head between leg-rest or beneath. ● Don't use the massage chair while wearing hair ornaments. ● Don't put hard things in pocket when take buttocks and huckle massage. ● Don't sit or press on the backrest when the chair is in reclining position. ● Don't massage on exposed skin or with thin clothes. It may cause skin ache. 	To avoid accidents or injury
	<ul style="list-style-type: none"> ● Keep the chair away from heat, wet, sharp things and corrosive, inflammable, explosive environment. ● To use electric blanket or other heating products on the chair is forbidden. 	To avoid fade, discolor on upholstery or hardened, deteriorated on PU.
	<ul style="list-style-type: none"> ● Keep the chair away from stove or other heating products. ● Don't put the chair directly under the sun or other high temperature places. 	
	<ul style="list-style-type: none"> ● Don't use the chair outdoor. 	

**Caution**


	<ul style="list-style-type: none"> ● Don't use or store the chair in the wet environment like nearing the swimming pool, or bathroom. ● Keep the chair and remote control dry all the time. ● Don't use the chair when your body is wet. 	To avoid leakage or electric shock.
	<ul style="list-style-type: none"> ● Make sure to install the earth wire in correct way (Incorrect install way may cause electric shock. If you don't know how to install, please find a professional electrician) ● Don't use plug adapter (If the enclosed plug doesn't match with socket, please don't change the plug by yourself and find a professional electrician to install a matched socket.) 	To avoid electric shock.
	<ul style="list-style-type: none"> ● Unplug the power before cleaning the chair. ● Unplug the power after using the chair. 	To avoid electric shock or fire hazard.
	<ul style="list-style-type: none"> ● Unplug the power when any damage on the power. ● Hold the plug to unplug the power but not try to pull the wire to unplug the power. 	
	<ul style="list-style-type: none"> ● Please turn off the power switch and unplug the power after massage. 	To avoid accidents
	<ul style="list-style-type: none"> ● Unplug the power immediately when electricity goes off. 	To avoid damage or accidents when power on again.

Grounding instruction

- ✧ ※The chair is equipped with a cord having an grounding conductor. The plug must be plugged into an appropriate socket that is properly installed and matched with the plug.

Replace fuse

- ✧ ※When the fuse is burned out, please turn off the switch and unplug the power. Opening the lid covered the fuse which is located in the switch box. Using a same size fuse to replace the old one, then cover the lid. Or you can find an electrician to replace the fuse.

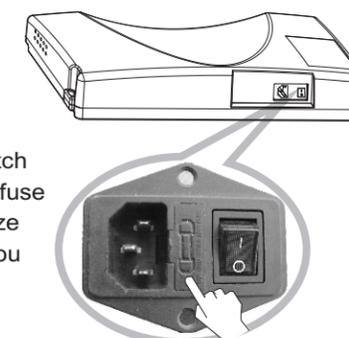

Earth wire is not allowed to be in these places

Gas pipe: probably cause explosion or fire.
Phone wire and lightning rod: probably cause electric shock and fire when thundering.
Water pipe: earth wire is useless where there are plastics

External structure




Preparation

Internal structure



Preparation

Explosion Diagram

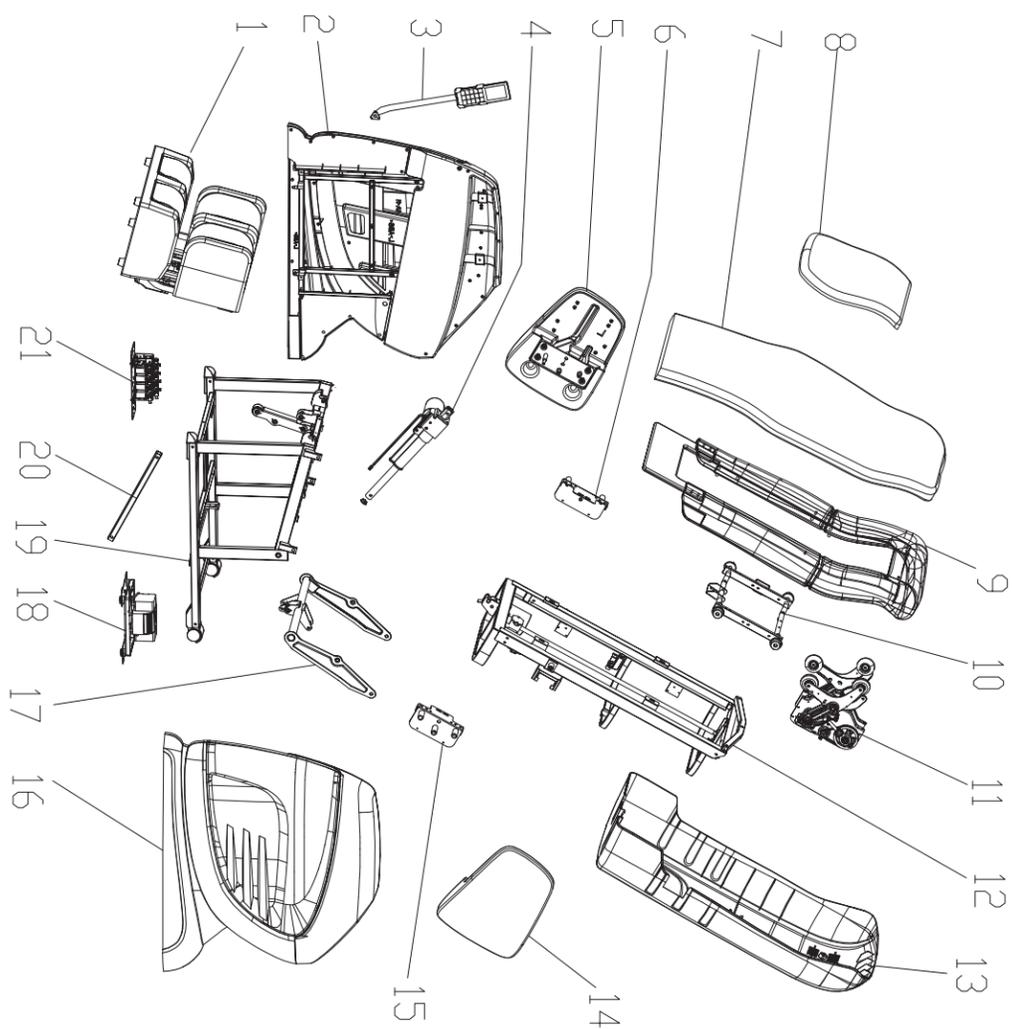

Preparation

Mechanical Components List

序号	名称	数量			
1	Legrest and footrest assembly	1			
2	Right armrest assembly	1			
3	Control bracket	1			
4	Electric actuator	1			
5	Right shoulder assembly	1			
6	Shoulder external fixation components (right)	1			
7	Back cushion	1			
8	Head pad	1			
9	Front backrest assembly	1			
10	Metal frame for massage machine	1			
11	Massage machine	1			
12	Backrest frame assembly	1			
13	Backrest cover	1			
14	Left shoulder assembly	1			
15	Shoulder external fixation components (left)	1			
16	Left armrest assembly	1			
17	Rotating frame	1			
18	Air pump assembly	1			
19	Seat frame	1			
20	Connecting rod	1			
21	Magnetic valve assembly	1			

Preparation

Installation (1)

Step 1. How to move the chair

Two or more people grasp left and right metal seat frame, take the chair body, back cushion, head pad, seat cushion controller bracket and electric power cable out of the packing box.


(as shown in the picture)



Main body of the chair

Hold the metal seat frame and take the chair out of the box.

Keep the foam piece in the original place between the backrest and seat while taking the chair out of the box in case PU leather gets broken due to moving rubbing.

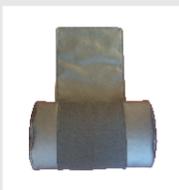

Head pad



Head pad



Seat cushion



Control bracket



Electric power cable



Back cushion



Right shoulder massage device




Left shoulder massage device

Step 2. How to assemble the backrest

Raise the backrest of the chair to 120 degree till "click sound" comes out which means the backrest is locked well with the buckle under the seat.



Note: When the backrest is leaning forward, open the seat cushion from here (please refer to the picture), pull up the lock for the backrest to lean forward, if you are not professional repairman, please don't operate this step.

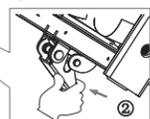

Installation (2)

Step 3 Assemble shoulder massage device

① Make the three axles on the shoulder massage device align to the three axle holes on the backrest. ②Unzip the shoulder massage device, insert the bolt screw tightly, zip it in the end. ③Repeat the above steps to install another shoulder massage device.


★ Disassemble shoulder massage device: Unzip the shoulder massage device, screwing out the screws; Turn the shoulder parts upwards, move the axle away from the axle hole, then disassemble shoulder massage device.



① Make the three axles on the shoulder massage device align to the three axle holes on the backrest, the longest axle on the shoulder device corresponds to the bottom axle hole on the backrest side.



②

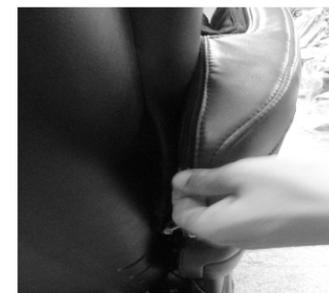

Installation (3)

Step4. Install the controller holder and put controller in.

Using three screw (M5*16) to fix the manual control to right upper arm.
Use the cover round holder cover the screws.
Then insert the manual control into the bracket.



Step5. Connect seat air pipe

Connect seat cushion with air pipe from the seat frame, then connect the velcro of seat cushion and seat frame.




Installation (4)

Step 6, Sorting out the zipper and snap fastener
Connect the back cushion, pillow with the backrest.



★Kind suggestion


Sitting on the chair, try to keep the bottom of pillow and ear in the same line, then connect it with back in this position. It will make body detection more accurate and provide better massage.

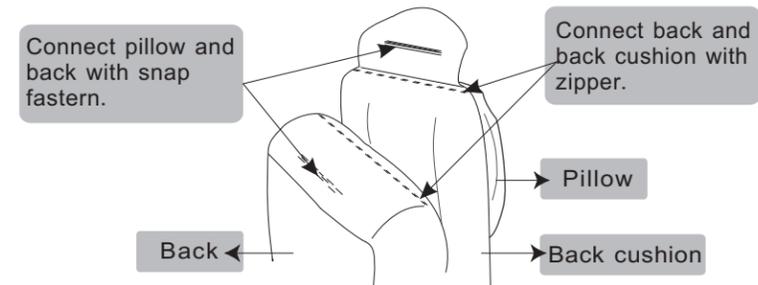

Step 7. After installation

After installation, please follow the steps below:

- 1 Insert the plug of remote control and power into the power box in bottom (as picture)
- 2 Insert the other end of power plug into the outlet.
- 3 Turn on the power switch (Power switch "O/I", "I" means on)
- 4 Press keys on remote control to choose massage program to enjoy the massage.

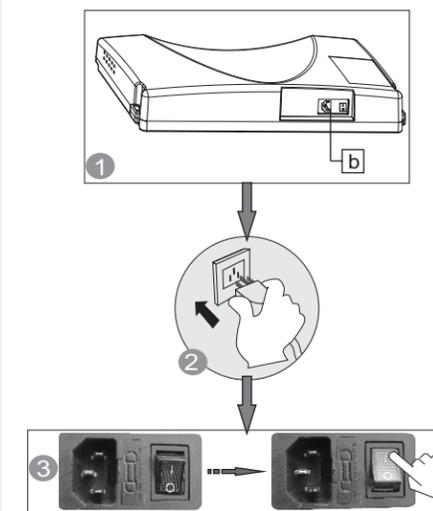

Environment



Caution



Don't use the chair in wet places like near swimming pool, bath room to avoid leakage or electric shock.



Don't use under sunlight or near stove or other heating place to avoid upholstery deteriorated.



Please use chair in flat place to avoid chair falling or noise or other unexpected trouble.

- Check places in front

Check the areas when back and leg extending, make sure no people, pets or other things around.



- Check power cable, plug

- (1) Clean the dust on plug on time. (using dry cloth).
- (2) No damage, re-do, bend, stretch, tie or press with heavy thing (especially can't be pressed by the chair) on power cable.



Check environment

- Space size for using the chair

- At least 40cm from wall.



- At least 40cm in front of the leg.

- Check surroundings

When the backrest or legrest reclines or extends, please pay attention to:

- (1) Check any person or things in every direction.
- (2) Check the space size.



How to move massage chair



Caution

- ★ Don't move chair when it's working in case of accidents or damage


How to move

- ※ Bind the cables and put them in the seat to avoid damage when moving.
- ※ If you plump the chair, it may damage the inner parts of the chair may get damaged.
- ※ Don't use the castor if the floor is wooden or other easy damaged material. It will need 2 or more persons to lift the chair.
- ※ Make sure the power is off and the plug doesn't touch the floor (It's better to let the chair in standing status to avoid the back crashing the floor)



Normally

- 1 Lift up the legrest to make the gravity fall on the castor. Then push chair to destination and set it down.

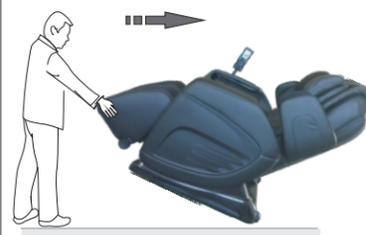

- ⚠ Please hold the footrest to move the chair as shown in the picture in case any scratches on the leatherette

- 2 Two person hold the two arms, lift and moving the chair to destination.



Adjustment before massage

- Before adjustment

- ★ Please make sure the massage head in right position before sitting on the chair (to avoid massage roller hitting your head)



- ※ Acupuncture point track testing confirm acupuncture point
If massage head can't detect when acupuncture point position detecting, shoulder position was default to the highest position.

- ★ Correct pose: Sit down and put the back closely to the backrest of the chair, user's head should lie down back to the head pad. During body height testing period, the massage roller will softly roll to the middle of the backrest before coming out to maximum, then go to the top to test shoulder height.



- ※ Avoid for the massage head detecting your acupuncture point position not right, otherwise you won't get good massage.

- Kind suggestion



- When you feel the intensity isn't strong enough, take off the pillow or cushion.

- Because the chair is heavy, it may ruin your floor in long time setting. Please put a blanket or other soft things to protect your floor.



Manual control



Massage operation

-  1. Press power button under standby status, massage chair enters ON status. Press Menu or any Auto massage program. In ON status, press power button, massage chair will turn off and return, massage machine stays in upper backrest, legrest stays in the lowest position, backrest stays in the highest position, controller will turn off after 2-3 seconds.
-  2. In ON status, press it to enter quick relaxation auto program, there is relative status screen display. Default massage time is 8 minutes. Press it again to stop.
-  3. Press this key to stop massage, press it once again to start massage
-  4. Press this button in Main interface to switch to menu interface, user can set different function for massage chair. Press it in Menu interface to return to main interface.
-  5. In On status, press it to switch an auto program each time ,and relative screen display is shown, 6 auto programs for choice, "Recovery" , "Pain relief" , "Stretch" , "Overall" , "Rocking relax" and "Helth care" .
-  6. In On status, press it to enter Air pressure program, press it to switch an auto program each time, and relative screen display is shown.
-  7. In ON status, press it to adjust air intensity with 3 grades, air intensity is adjusted by each press, and relative screen display is shown.
-  8. Direction button has 4 buttons,they are Up,Down,Left and Right.In Menu interface it is navigator button, user can press Up or Down to choose menu fast, and press Left or Right to enter previous or next menu, What' s more, Up and Down button can be used to adjust massage position in acupuncture point tracking detecting finish and point massage.
-  9. In ON status and under Menu mode, press it to confirm menu choice and enter next menu.
-  10. Keep pressing this button, backrest go up, legrest go down, release this button, backrest and legrest stop moving. When adjusting sitting position, it shows seat position adjusting... in controller screen.
-  11. Keep pressing this button, backrest go down, legrest go up, release this button, backrest and legrest stop moving. When adjusting sitting position, it shows seat position adjusting... in controller screen.
-  12. In On status, press this button, massage chair will go into zero-gravity position, press it again to return to default position.

Massage operation

LCD display graphic description

Main interface



Press "on/off" button when power on, a buzzer beeps once short and brittle. LCD screen display boot screen, after 2-3 second it display the main interface and enter standby mode. Main interface show massage position, methods, time, current massage status and so on, all information will be updated timely when the chair start working.

- Massage function
- Kneading
 - Tapping
 - Knocking
 - Pressing
 - Zero gravity
 - Massage time left




- ▬▬▬ Massage speed (slow, mid, fast)
- ▬▬▬ Massage width (narrow, mid, wide)
- ▬▬▬ Air intensity (weak, mid, strong)

Menu interface



Press the "menu" button in the main interface, it will switch to the menu interface, user can set all massage function. Press the "menu" button in the menu interface will exit and return to the main menu interface.

Auto : press "menu" button after the cursor is the default in automatic massage program menu item, press the "OK" button, there are five auto massage modes can be chosen.

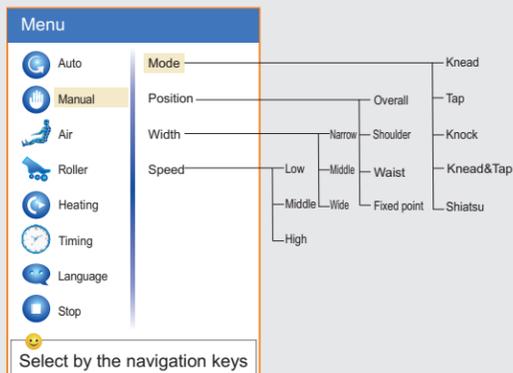


Manual : press "menu" key, select manual massage by the navigation key, press the "OK" button, there are four options: mode, position, width adjustment, speed adjustment : select the corresponding function, then press the "OK" button to enter the next menu, select the desired massage method, position, width and speed regulating function.

When massage chair in manual fixed point work modes, users can use up/down button to adjust massage position, information for "massage position adjusting..." will be show at the bottom of the display screen.

Massage operation

LCD display graphic description

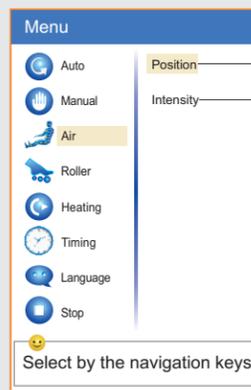

Air : After press MENU button, choose "air pressure" through navigation key, Press "OK" button , there are two kinds of air pressure related function: Position and Intensity; Choose the corresponding function, then press "OK" to enter next menu and choose required function.

Roller : press "menu" button, select foot roller massage by the navigation key, press the "OK" button, Press "OK" button, choose rolling massage mode.


Heat : in menu mode, use navigator key to look for waist heat, press on or off to choose heat on or heat off

Timing : press "menu" key, through the navigation key to select menu item timing , can choose massage time : 30 minutes, 20 minutes, 10minutes, the default time is 20 minutes, the chair will turn off automatically after time is up.

Language : press "menu" key, select language through the navigation key, then select the needed language.


Stop : press "menu" key, select stop through the navigation key, it will stop massage, and back to main interface.

After Power is on, the first time to use auto program, massage chair will start acupuncture point position detecting, and "acupuncture point tracking ..." will be shown on screen display. After acupuncture point tracking detecting is finished, buzzer will tick 5 times (once every second), meanwhile, the massage machine is located in shoulders position, if user think the position is not proper, press Up and Down to adjust, when adjusting, "acupuncture point adjusting..." is showing on screen display. Default working time is 20 minutes. (Massage chair will turn off and return after the massage is over)

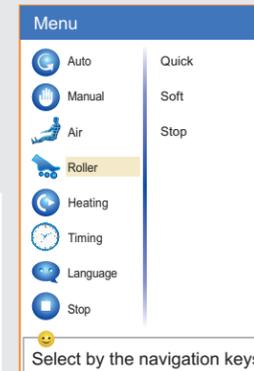

Massage operation

ON/OFF

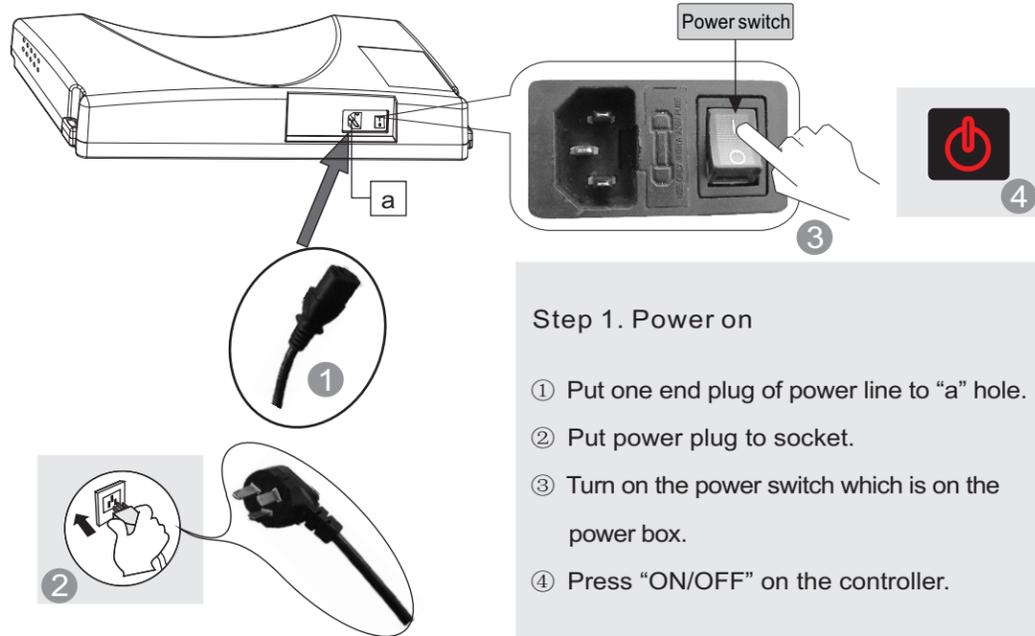

Step 1. Power on

- ① Put one end plug of power line to "a" hole.
- ② Put power plug to socket.
- ③ Turn on the power switch which is on the power box.
- ④ Press "ON/OFF" on the controller.

Step 2: Massage

- ① Operate according to controller operation instruction.
- ② Time is up, shut down massage.



Step 3: Turn off the power

- ① When the chair is ON, press "ON/OFF" button to turn off the massage chair and it will restore automatically.
- ② Turn off the power switch on the bottom of the back of chair after the rollers stay at storage position.
- ③ Unplug the power plug from the socket.



Warning

- ① Before using, pull up the pillow, back cushion etc. to check if the PU material is broken, if there is damage, Please stop using, unplug the chair, and ask for repairment.
- ② Only plug to proper socket.



Attention

- ① Before sitting down, make sure there is no things pinched between armrest, legrest, backrest and seat of the chair.
- ② Make sure the chair is in returned position before using.
- ③ Do not stand on it when using the massage chair.

Function Illustration

● Backrest up and legrest down

Press this button on controller, backrest go up and legrest go down. Release this button, backrest and legrest stop moving, the angle is locked in this position.



● Legrest extend and retract

Legrest can extend or retract within a range of 20CM.



● Legrest up and backrest down

Press this button on controller, backrest go down and legrest go up. Release this button, backrest and legrest stop moving, the angle is locked in this position.



● Arm massage

The armrest has pulse pressure massage function, put arms between arm airbags before the airbags filled with air. It will perform pinch massage technique to arms and hands after airbags inflating.



● Shoulder massage

Multi-layer extrusion airbags is deployed on both sides of the shoulder. Achieve pressure shoulder stretch massage, relieve shoulder fatigue.



● Zero gravity

Press the button "Zero", backrest go down and legrest raise up, the chair go into Zero gravity status automatically. Press the Button again, the backrest and legrest stop moving, and it will lock at this position. However to change the angles, we can keep comfortable massage status all the time.



⚠ Warning

When press the angle adjusting button, don't put hand or arm in the gap between shoulder and armrest, in order to avoid damage.



1. Please make sure there is no people, no pets or any stuff within range for backrest to avoid unexpected hurt or property loss.
2. When adjusting the backrest, be careful the gap between the backrest and armrest.
3. When the legrest is moving, please do not stand on or leave from the seat, do not stop the legrest from moving in case of causing damage to the chair.
4. Please refer to control instruction for more details.

Restoration of massage chair

- 1 Make sure there is no people or pet or anything within the range of backrest and legrest.



- 2 Make backrest and legrest restore to the normal position by controller.
(Or press Backrest adjust button to make backrest and legrest restore to the normal position)




- 3 After restoring, put the controller back into the controller holder.




- 4 Put the back cushion and pillow pad on to the right position.

- 5 Switch off the power.

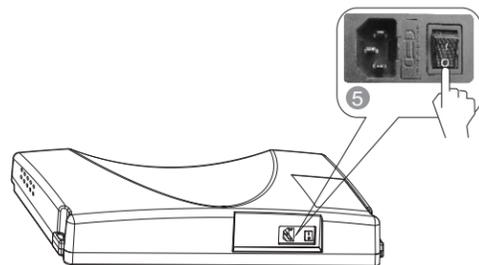

- 6 Unplug the electricity from the socket.
(Attention: hold the right plug position)

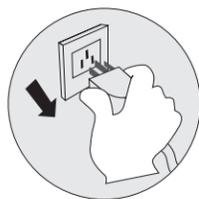

- 7 Clean the chair if necessary.



Attention!

- After using, press ON/OFF to end all the massage method and wait for restoration. Make sure to cut the electricity to avoid any unexpected hurt to children or pets.
- After using, turn off the power, unplug the electricity. Avoid moisture, dirt and aging of insulation in case of causing damage.

Synthetic leather

● General cleaning

Clean with tidy, soft and dry cloth slightly.




(Chemical or medical cleaning products are forbidden)



● If leather is dirty

- Dip soft cloth into 3-5% neutral detergent solution, squeeze it then clean the stain on the surface gently.
- Wash cloth, squeeze it and wipe with detergent.

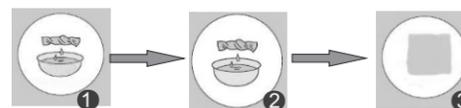

- Clean it softly with dry cloth.
- Allow to dry (air blower is prohibited)

Cloth material

- Dip cloth into weak neutral detergent, squeeze it and clean the material surface.

- (Diluent, gas, alcohol are prohibited)
- Brush the dirty parts with weak neutral detergent. (avoid excessive brushing)
Excessive cleaning may cause damage to cloth material.

- Wipe with water-dipped but squeezed cloth.
- Allow to dry.



Plastic parts

- Dip cloth into neutral detergent, squeeze it and clean.
- Wipe with water-dipped but squeezed cloth.
- Allow to dry.



Clean the controller

- Wipe with soft cloth after squeezing.
- Allow to dry.

Propellant, diluent or alcohol are prohibited.



About cleaning

- Before cleaning, unplug the electricity, don't touch plug with wet hands.
- Make sure hands are dry when plugging to the socket.
- If failed to follow the mentioned instruction, user may get electric shock or burned.

FAQ and Troubleshooting

Trouble	Possible reason	Troubleshooting
The massage chair can't work after power is on.	Check if the control is on or not.	Open control switch
	Not choose massage function	Choose massage function
	Check whether the plug and socket in poor contact.	Make sure plug and socket are connected.
	Wire or plug has damaged.	Please ask manufacturer or appointed or authorized maintenance man for checking and replacement.
	The fuse burned out.	Replace fuse of same specification.
	Something wrong with inner circuit	Please ask appointed or authorized maintenance man for checking and repairing.
There is beep sound inside the massager while working.	It is sound of air pump, motor and mechanical parts.	Normal sound, no need to pay attention to it.
Different height of the left roller and right roller.	Massage rollers work on shifts.	Normal, no need to pay attention to it.
Sound suddenly becomes much louder while working.	Check if it is due to long time working.	Turn off the power to have the massager rest for over half an hour before next turn working.
	Inner parts fretted away due to long time work. It's normal	Please ask appointed or authorized maintenance man for checking and repairing.
Massage function suddenly stops while working.	Suddenly power-down	Connect the electricity power again
	Setting time is up.	Turn off the power to have the massager rest for over half an hour before next turn working.
Neither the backrest nor the legrest can be lay down.	Check if there are any obstacles; the chair overloads; for the sake of safety, sometimes the massage machine may automatically stop working.	Remove the obstacles; turn off the power to have the massager rest for over half an hour before next turn working.
The chair can't return to its original position.	Excessive operation; setting time is up.	Please ask manufacturer or appointed or authorized maintenance man for checking and replacement.
Electric wire or plug give out heat	It may work overload	Turn off the power to have the massager rest for over half an hour before next round working.
	Other abnormal case	Stop using, and ask appointed maintenance man for help.



Note: If the troubles are not caused by above reasons, please turn off the power switch and cut off the power, then contact the manufacturer or service agent to repair your chair.

Technical data

Name	Massage Chair	Power	180
Model	SERENITY ZA-151	NW/GW	76/88kgs
Input voltage	<input type="checkbox"/> AC120V, 60Hz <input type="checkbox"/> AC220V, 50/60Hz <input type="checkbox"/> AC240V, 50Hz	Packing size	1260*780*810mm (L×W×H)
Safety design	Class I		